



Electrochemical Analytic Chip

Electrochemical analytic chip, version 6, encapsulated on a ceramic board for electrochemical measurements in liquids.

Service

The measurement of voltage and current is essential in electrochemical analysis for determining chemical compounds in solution and their concentrations. Electrodes made of gold or platinum, in combination with reference electrodes such as Ag/AgCl, are often used for this application. A typical disadvantage of conventional measurement setups is the need for at least 20 ml of analyte, which requires large electrodes. As a consequence, the electrodes can be expensive and have to be cleaned periodically with substantial effort. Additionally, contaminated electrodes might cause measurement errors, where the replacement of electrodes made of precious metals is costly and environmentally unfriendly. The analytic chips developed by IPMS can overcome these issues and are especially suited for potentiometry and voltammetry measurements. The chips, with an area of 5 x 5 mm², offer different electrode geometries, sizes, and materials for the working, counter, and reference electrodes. Precise electrode geometries and dimensions provide another advantage for electrochemical measurements. The chips with these electrodes are easy to handle and require only a very small amount of precious metal due to their resource-efficient fabrication. Analysis with the chips can be performed directly with


a measuring adapter, or they can be bonded to a ceramic board. The size of the ceramic boards is 30 x 9.3 x 0.63 mm³ and the boards have six contact pads, allowing for the use of up to two individual 3-electrode measuring cells on one chip. A sample size of one droplet (~20 µl) containing the analyte is sufficient. The setup is compatible with aqueous solutions (pH 2... pH 10), as well as with most organic solvents. Available electrode materials are platinum, gold, silver and copper, with thicknesses from 1 µm to 3 µm. Customized chips, sensors, electrode geometries, and electrode materials are possible, thereby covering a broad spectrum of applications.

Contact


Dr. Olaf Rüdiger Hild
Tel. 49 351 8823-450
olaf.hild@ipms.fraunhofer.de

Fraunhofer Institute for
Photonic Microsystems IPMS
Maria-Reiche-Str. 2
01109 Dresden
Germany

www.ipms.fraunhofer.de

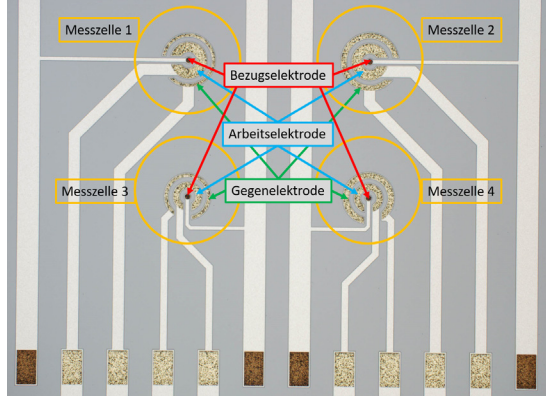


Close up of chip no. 6 with 4 different 3-electrode set ups

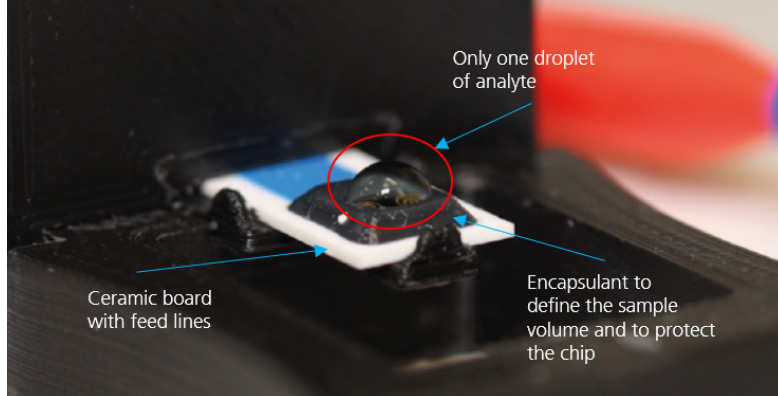


Chip module with one droplet of analyte in a measuring adapter.


Available chip geometries and electrode sizes

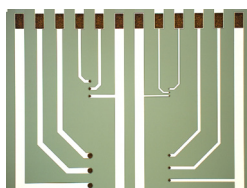


Electrode arrangement 1	Electrode geometry	Electrode distance
Messzelle 5 – 3x_D20 μm _50 μm	20 μm , rund	50 μm
Messzelle 6 – 3x_D50 μm _100 μm	50 μm , rund	100 μm
Messzelle 7 – 3x_D50 μm _200 μm	50 μm , rund	200 μm
Messzelle 8 – 3x_D100 μm _200 μm	100 μm , rund	200 μm

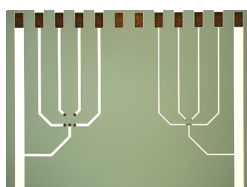


Electrode arrangement 2	Electrode geometry	Electrode distance
Messzelle 13 – 5x_D20 μm _80 μm	20 μm , rund	80 μm
Messzelle 14 – 5x_D50 μm _160 μm	50 μm , rund	160 μm

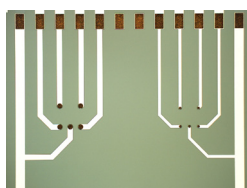


Electrode arrangement 3	Electrode geometry	Electrode distance
Messzelle 15 – 5x_D50 μm _320 μm	50 μm , rund	320 μm
Messzelle 16 – 5x_D100 μm _320 μm	100 μm , rund	320 μm

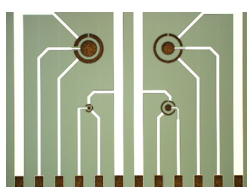

Electrode arrangement 4	Electrode geometry	Electrode distance	Real size μm^2
Messzelle 17 – 3x_D60 μm _30 μm	60 μm , rund	30 μm	2800, 4300
Messzelle 18 – 3x_D130 μm _60 μm	130 μm , rund	60 μm	13000, 13000
Messzelle 19 – 3x_D330 μm _60 μm	330 μm , rund	60 μm	85500, 85500
Messzelle 20 – 3x_D260 μm _120 μm	260 μm , rund	120 μm	53000, 53000

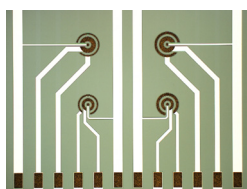

Electrode arrangement 5	ID, d	Radius	Real size μm^2
Messzelle 25 – WE-CE-RE	50 μm , 50 μm	75 μm , 125 μm , 175 μm , 202 μm	27000, 27000
Messzelle 26 – WE-CE-RE	50 μm , 50 μm	75 μm , 125 μm , 175 μm , 223 μm	27000, 54000
Messzelle 27 – WE-CE-RE	50 μm , 50 μm	75 μm , 159 μm , 209 μm , 232 μm	54000, 27000
Messzelle 28 – WE-CE-RE	50 μm , 50 μm	75 μm , 159 μm , 209 μm , 253 μm	54000, 54000

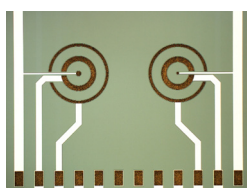

Electrode arrangement 6	ID, d	Radius	Real size μm^2
Messzelle 29 – WE-CE-RE	100 μm , 200 μm	250 μm , 350 μm , 550 μm , 603 μm	180000, 180000
Messzelle 30 – WE-CE-RE	50 μm , 200 μm	225 μm , 333 μm , 533 μm , 588 μm	180000, 180000

*WE = Working electrode, CE = Counter electrode; RE = Reference electrode
 Parts of the underlying research work were funded by EFRE-Funds of the European Union
 and by the SAB (Sächsische Aufbaubank) under grant agreement number: 3784.